



GB119



GB149

Owner's Manual

Buzzaround Luxury™



Thank you for purchasing your Buzzaround Luxury™ Scooter. We are honored you have chosen Golden Technologies for your mobility needs.

Buzzaround Luxury™
Models: GB119 and GB149

Standard Packing List:

- 1 Owner's Manual
- 2 Battery Packs
- 1 Front Basket
- 1 Seat with Armrests
- 1 Safety Flag
- 1 Interchangeable Color Panel Set
- 2 Keys
- 1 Battery Charger



DISCLAIMER

We ask that you read this manual completely before operating your new Buzzaround Luxury™. Golden Technologies, Inc. is not and cannot be held responsible for any damage or injury incurred due to improper or unsafe use of the Golden Buzzaround Luxury™ scooter. Golden Technologies, Inc. specifically disclaims responsibility for any bodily injury or property damage that may occur during any use that does not comply with applicable federal, state, or local laws or ordinances.

Intended Use: The intended use of the Buzzaround is to provide mobility to persons limited to a seated position that are capable of operating an electric scooter.

CONTENTS

I. INTRODUCTION	5
Features and Benefits	6
Accessories for the Buzzaround Luxury™ Scooter	7
II. SAFETY	8-12
III. EMI/RFI	13
IV. YOUR GOLDEN Buzzaround Luxury™	14
Specifications	15-16
V. OPERATION	17
Control Panel	17-18
Control Options	18
USB, Fuse Panel, and Seating	19-20
Freewheel Mode	21
VI. RIDING YOUR Buzzaround Luxury™ Scooter	22
Mounting	22
Basic Driving	23
Steering in Tight Spots	24-25
VII. DISASSEMBLY	26-27
VIII. ASSEMBLY	28-31
IX. TRANSPORTING YOUR Buzzaround Luxury™	32
X. CARE AND MAINTENANCE	33
Routine Maintenance	33
Cleaning Your Buzzaround Luxury™ Scooter	33
Batteries and Charging	34
If Your Buzzaround Luxury™ Does Not Operate	35
Proper Tire Inflation for (Optional Pneumatic Tires)	36
Diagnostics	36-38
XI. WARRANTY	39-40
Warranty/Registration Card	41

I. INTRODUCTION

Congratulations on the purchase of your new Buzzaround Luxury™ scooter. The Buzzaround Luxury™ combines cutting edge technology with attractive designs that are also highly functional in today's world. We at Golden Technologies, Inc. know that you have chosen a scooter that will give you years of dependable operation and also will enhance the quality of your life by providing you with the mobility to experience an active daily lifestyle.

Even though your new Buzzaround Luxury™ is both user-friendly and designed for maximum maneuverability in even the tightest spaces, we ask that you please read, understand and follow all of the instructions and suggestions in this manual before you operate your scooter for the first time. The safe use of your new scooter is very important to us.

If you feel that you do not understand the instructions and suggestions presented in this owner's manual, or if, for any reason, you do not feel capable of performing the activities necessary to assemble, disassemble, operate, or maintain your Buzzaround Luxury™, please contact your local Golden Technologies, Inc. provider or call Golden Technologies, Inc. Technical Support Services at (800) 624-6374.

Golden Technologies, Inc. cannot be held responsible for personal injury or property damage resulting from the unsafe or the improper use of any of our broad range of health and personal mobility products. Also, Golden Technologies, Inc. cannot be held responsible for personal injury or property damage resulting from attempts to follow instructions, suggestions, and guidelines presented in this owner's manual.

Our Research and Development Department, our Quality Control Department, and our Engineering Department have used the latest product specifications and the latest product design information to manufacture your Buzzaround Luxury™. Golden Technologies, Inc. reserves the right to implement changes into our product lines when those changes become desirable or necessary. If changes are implemented into our product line, there may be minor differences between the product you purchased and the illustrations and instructions in this owner's manual.

Please fill out and mail the enclosed warranty registration card. Registration is also available @ www.goldentech.com. We at Golden Technologies, Inc. would appreciate hearing about the dependability of your Buzzaround Luxury™ and about the convenience of mobility it has provided for you. We would also appreciate hearing about the service you received from your local Golden Technologies, Inc. provider or representative.

Golden Technologies, Inc.
401 Bridge Street
Old Forge, PA 18518

Phone: (570) 451-7477
Fax: (570) 451-7494
Toll free: (800) 624-6374

For more information about our products and services or to send us your comments, please visit our website at www.goldentech.com

I. INTRODUCTION

FEATURES AND BENEFITS

Your Buzzaround Luxury™ scooter has been designed with your comfort in mind. It has also been designed to be a safe method of transportation when safety and operating instructions are followed. Please review the features and benefits of certain components of the Buzzaround Luxury™ scooter below.

Delta Tiller Handle:

Offers rubber coated hand rests for comfort and dual throttle control levers for easy right or left hand use.

Headlight and Taillight On/Off Switch:

The switch is located in the center of the tiller under the control panel.

Basket:

The Buzzaround Luxury comes standard with a removable storage basket.

Tiller:

Easy to use tiller offers multiple adjustments to find your most comfortable driving position. The tiller angle adjustment lever is located on the right side of the tiller, and locks automatically in place when released. A battery charging port, USB charging port, and fuse panel are conveniently located in the center of the tiller under the control panel.

LED Headlight:

The Buzzaround Luxury comes standard with a LED angle-adjustable headlight and LED taillight, which is designed to make it easier to see and be seen.

Control Panel:

Easy to read display features from top to bottom: the LED battery gauge; oversized speed control knob; key switch to turn the unit on and off; and the red horn button.

Adjustable Armrests:

Your Buzzaround Luxury scooter features width and angle adjustable armrests to find the most comfortable position as well as to make transfers easier.



Seat Swivel:

A standard feature, the seat can swivel to the left or right for easy transfers or to sit more comfortably at a table or desk.



Easy Access Battery Packs:

Located under the seat, the battery packs have a wireless connector system for easy disassembly and reassembly.

Interchangeable Color Panels:

Located on the front and rear fenders of the scooter. Standard red and blue panels can be quickly changed by simply pushing the panels up from the access hole located underneath the plastic shroud.

Figure 1. Buzzaround Luxury™ Features (Model: GB119)

I. INTRODUCTION

ACCESSORIES FOR THE BUZZAROUND LUXURY™ SCOOTER:

A variety of accessories are available for your Golden Buzzaround Luxury™ Scooter. Please contact your Authorized Golden Technologies Provider for more information or to order.

- Cane Holder - Tube
- Quad Cane Holder
- Walker Holder
- Oxygen Tank Holder
- Cup Holder
- Deluxe Pack N' Go
- Forearm Crutch Holder
- Mirror
- Safety Flag
- Scooter Cover
- Cell Phone Holder

*Some accessories require a mounting bracket or mounting clips. Please check with your provider to find out if the accessory you would like to order requires a mounting bracket or mounting clips.

Safety Guidelines

The symbols below are used in the owner's manual and/or on the product to identify warnings and to notify you of important safety information. Make sure you read and understand them completely.



WARNING! Indicate(s) a hazardous situation which, if not avoided, could result in death or serious injury. This icon is represented as a black symbol on a yellow triangle with a black border.



MANDATORY! Failure to perform mandatory actions can cause personal injury and/or equipment damage. This icon is represented as a white symbol on a blue dot with a white border.



PROHIBITED! Indicate(s) an action(s) that should not be performed under any circumstance. Failure to adhere can cause personal injury and/or equipment damage. This icon is represented as a black symbol with a red circle and red slash.

Safety Symbols

The symbols below identify warnings, mandatory actions, and prohibited actions. They can be found in the owner's manual and/or on the product. Make sure you understand all safety labels before operating the product. Do not remove safety labels from the product.



Read and follow the owner's manual!



Explosive conditions exist!



Indoor use only!



Use only AGM or Gel-Cell batteries!



Avoid transmitters!



Hot surface!



Pinch/Crush Hazard!



Do not place into trash!



Electrical hazard!



Corrosive chemicals inside battery!



Recycle



Do not use personal electronic devices!

II. SAFETY

Your Buzzaround Luxury™ is a battery-operated personal mobility vehicle. Please exercise caution and consideration when you are operating it. Driving your Buzzaround Luxury™ carefully and thoughtfully will help ensure your personal safety and the safety of other people.

NOTICE

Before learning to operate your Buzzaround Luxury, have your authorized Golden Technologies, Inc. provider determine if it is advisable for you to practice getting on and off your scooter and operating it in the presence of an attendant.



MANDATORY! Read completely and understand this owner's manual before assembling, operating, transporting, or disassembling your Buzzaround Luxury scooter.

BEFORE GETTING ON YOUR BUZZAROUND LUXURY

- Check to be certain that the power is turned off. See section V. "Operation" in this manual. This will eliminate the possibility of accidentally activating the paddle controls and causing injury to yourself and/or others.
- Check to be certain that your Buzzaround Luxury's freewheel lever is in the engaged position. See section V. "Operation" in this manual.
- Flip up the armrest(s).



WARNING

When getting on or off your Buzzaround Luxury, keep your weight toward the middle of the deck. Putting most or all of your weight on the edge of the deck may cause an unstable condition.

GETTING ON YOUR BUZZAROUND LUXURY

- Position the seat for safe and easy mounting. Ensure that the seat is facing forward. See Section V. "Operation."
- Return tiller to full upright position if necessary.
- Carefully place one foot on the approximate center of the deck and seat yourself comfortably and securely on the seat.
- Flip down the armrest(s).

GETTING OFF YOUR BUZZAROUND LUXURY

- Make certain that the power is turned off and the key is removed.
- Return tiller to full upright position.
- Flip up the armrests.
- Carefully place one foot on the ground, transfer your weight to that leg, and slowly come to a standing position.
- Step away from the scooter.

MAXIMUM WEIGHT

Your Buzzaround Luxury has been rated to a maximum weight capacity (operator and anything else being carried onboard) of **375 pounds**. Exceeding the maximum weight rating will void the warranty.



WARNING

Exceeding the maximum weight rating may result in injury to yourself and/or others.

II. SAFETY

DRIVING ON INCLINES

- For maximum stability, lean forward in your Buzzaround Luxury's seat while proceeding up ramps, inclines, curbs, or any low rise. See figure 2.
- Drive with caution when attempting to negotiate any incline, even handicap access ramps.
- **Always** climb or descend an incline by driving straight up or straight down the face of the slope. See figure 2.
- **Do not** traverse or drive across the face of an incline in any direction. See figure 3.
- **Do not** attempt to negotiate an incline that is covered with snow, ice, salt, cut or wet grass, leaves, or any other potentially hazardous material.
- **Do not** drive your scooter in reverse down an incline.
- Try to keep your Buzzaround Luxury moving when climbing an incline. If you do come to a stop, restart and accelerate slowly and carefully.
- **Do not** try to descend or climb a slope whose gradient is greater than recommended.

*** Buzzaround Luxury maximum recommended incline is 6°.**



Figure 2. Going Up an Incline



Figure 3. Traversing an Incline



WARNING

If, while you are driving down a slope, your scooter starts to move faster than you feel is safe, release the throttle control lever and allow your Buzzaround Luxury to come to a stop. When you feel that you again have control of your scooter, push the throttle control lever forward and continue safely down the remainder of the slope.

DRIVING DOWN A DECLINE

- Lower speed setting.
- Whenever it is safely possible, drive forward down any ramp, low rise, or incline.

Golden Technologies, Inc. does not recommend that you drive your scooter in reverse down any incline, ramp, curb, or low rise. Backing down any slope may create a very hazardous situation. If it is necessary to back down a slope follow one of two procedures on page 10.

II. SAFETY

CAUTION

When using either of the following procedures to back down a slope, Golden Technologies, Inc. strongly advises that you have the assistance of another person.

Procedure 1 - Under Power

1. Set the key switch to the “Off” position. See section V. “Operation.”
2. Turn the speed control dial fully counterclockwise to its slowest setting.
3. **Dismount your Buzzaround Luxury™.**
4. Set the key switch to the “On” position. See section V. “Operation.”
5. While standing alongside your scooter, carefully operate the controls at the lowest speed setting. See section V. “Operation.”
6. Slowly and carefully back your scooter down the incline.
7. Carefully remount your Buzzaround Luxury and resume normal operation.



WARNING

When your scooter is in freewheel mode, the park brake is released. If the batteries are dead, the motor braking system will not work. The weight of the scooter on a slope may cause you to lose control of the scooter. If you do not feel capable of manually controlling your scooter down the slope, request assistance or do not attempt to use this procedure.

Procedure 2 - Emergency (No Power)

1. Set the key switch to the “Off” position. See section V. “Operation.”
2. **Dismount your Buzzaround Luxury™.**
3. Set the freewheel lever to “Disengaged.”
4. While standing alongside the scooter, carefully manually maneuver it down the slope.
5. When you have reached a level surface at the bottom of the slope, position the freewheel lever in the engaged position. See section V. “Operation.”
6. Carefully remount your Buzzaround Luxury, and set the key switch to the “On” position to restore power.

Motor Braking System

Your scooter is equipped with a system that uses the motor to aid in braking. This motor brake system is designed to work when the key is in either the on position or the off position.

When the key switch is in the on position, the freewheel is in the engaged position, and the scooter is under power, the motor will help slow down the scooter as soon as you take your hand off of the throttle lever.

When the key switch is in the off position and the freewheel lever is in the disengaged position, the motor brake system will keep you from pushing the scooter too fast (i.e. down a slope).

You may notice this when pushing the scooter. The scooter will move freely until you reach a certain speed. You will then encounter some resistance as the motor brake system is activated.

II. SAFETY

MEDICATION

Always check with your physician to determine if any of the medications you are taking may affect your judgment and/or your ability to operate your Buzzaround Luxury. Also check with your physician concerning your physical ability to operate a scooter.



WARNING

Do not connect or allow anyone except an authorized Golden Technologies, Inc. provider to connect any electrical or mechanical device to your Buzzaround Luxury. Unauthorized accessories will void the warranty and may cause injury.

NOTICE

Please remember that while on your Golden Buzzaround Luxury, you are a motorized pedestrian. You must observe and obey all pedestrian rules and regulations for the locale in which you are riding.

Please use your Buzzaround Luxury™ scooter often and allow it to expand the horizons of your daily life. The more mobility your Golden Buzzaround Luxury brings to you, the happier you will be, and the more satisfied we at Golden Technologies, Inc. will be! But as with all things, especially motor operated vehicles, observing a few rules and safety considerations will help ensure safe scooter operation. So please follow the rules below.

RULES FOR USE AND OTHER SAFETY CONSIDERATIONS



WARNING

Read completely and understand this owner's manual before assembling, operating, transporting, or disassembling your Golden Buzzaround Luxury.

- **Always** operate your Golden Buzzaround Luxury with thought, care, and safety.
- **Do not** attempt to use your Buzzaround Luxury on an escalator. Always use an elevator.
- **Do not** carry passengers under any circumstances.
- **Do not** mount or dismount your Golden Buzzaround Luxury unless the brake is engaged.
- **Always** make sure that the key switch is set to "Off" before mounting or dismounting your scooter.
- **Do not** back your Golden Buzzaround Luxury down an incline or across an uneven surface while seated.
- **Do not** turn your scooter suddenly at full speed.
- **Always** make sure the seat is locked in the forward facing position before operating your Buzzaround Luxury.
- **Always** come to a full stop before changing direction from forward to reverse or from reverse to forward.
- **Do not** operate your Buzzaround Luxury where you could not safely or legally walk.
- **Do not** climb ramps or curbs that exceed your scooter's capacity.
- **Always** approach curbs and inclines straight on.
- **Always** be aware of and careful near mechanical pinch points especially when assembling and disassembling your scooter.

II. SAFETY

RULES FOR USE AND OTHER SAFETY CONSIDERATIONS - continued



WARNING

- **Never** sit on your Buzzaround Luxury when it is being transported.
- **Always** fasten down your Buzzaround Luxury securely with an approved tie-down system while transporting your scooter.
- **Never** operate your scooter if it is not functioning properly.
- **Always** use caution when driving on soft or uneven surfaces such as grass or gravel. Also use caution on decks where there is no railing.
- **Never** drive on the roadway, except when you must cross the street.
- **Always** cross streets at intersections and use crosswalks or the most direct route, making sure that your path is clear and that you are visible to motor traffic.
- **Never** drive your Buzzaround Luxury up or down a step or curb that is higher than the ground clearance listed on the specifications chart (see pages 15-16).
- **Never** back up or down a step or curb.
- **Do not** drive your Buzzaround Luxury in icy or salted conditions.
- **Never** operate your Buzzaround Luxury while you are under the influence of alcohol.
- **Do not** operate or store your Buzzaround Luxury where it will be exposed to rain, snow, mist, and below freezing temperatures.



WARNING

Should you encounter a problem with the operation of your mobility device such as a mechanical or electrical problem immediately stop using the product. Disconnect the batteries then contact the provider from whom you purchased your mobility device for evaluation of the problem and repair if needed.

III. EMI/RFI

The rapid development of electronics, especially in the area of communications, has saturated our environment with electromagnetic (radio) waves that are emitted by television transmitters, cellular phones, citizen's band radios (CBs), amateur radios (ham radios), wireless computer links, microwave transmitters, paging transmitters, etc. These electromagnetic (EM) waves are invisible and increase in strength the closer one gets to the source of transmission. When these energy waves act upon electrical devices and cause them to malfunction or to function in an erratic or uncontrolled manner, they are referred to as Electromagnetic Interference (EMI) or Radio Frequency Interference (RFI).

EMI/RFI AND YOUR BUZZAROUND LUXURY™

All electrically powered vehicles, including scooters are susceptible to Electromagnetic Interference/Radio Frequency Interference (EMI/RFI). This interference could result in abnormal, unintended movement of your Golden Buzzaround Luxury™.



WARNING

Unintended movement or brake release could cause an accident or injury.

The FDA has determined that each make and model of scooter can resist EMI/RFI to a certain level. The higher the level of resistance, the greater the degree of protection from EMI/RFI—measured in volts per meter (V/m). The FDA has also determined that current technology is capable of providing 20 V/m of resistance to EMI/RFI, which would provide useful protection against common sources of interference. This product has been tested and has passed an immunity level of 20 V/m.

EMI/RFI RECOMMENDATIONS



PROHIBITED! Do not turn on or use hand-held personal electronic communication devices such as cellular phones, walkie-talkies, and CB radios while your scooter is turned on.

- Be aware of any nearby transmitters (radio, television, microwave, etc.) on your intended route and avoid operating your scooter close to any of those transmitters.
- Turn off the power if your Buzzaround Luxury is going to be in a stationary position for any length of time.
- Be aware that adding accessories or components or modifying your scooter may make it more susceptible to EMI/RFI.
- If unintended movement or park brake release occurs, turn your scooter off as soon as it is safe to do so.
- Report all incidents of unintended movement or park brake failure to Golden Technologies, Inc. Technical Service at (800) 624-6374.



WARNING

Turn off your scooter as soon as it is safely possible if unintended or uncontrollable motion occurs or if unintended park brake release occurs.

IV. YOUR GOLDEN Buzzaround Luxury™



Tiller Controls:

Tiller angle adjustment lever, headlight/taillight switch, XLR battery charging port, USB charging port, and fuse panel.



Figure 4. Your Golden Buzzaround Luxury™ (Model GB149)

FOR YOUR RECORDS

Please fill in your Buzzaround Luxury's information below. This information will be useful in the event that you ever need to contact Golden Technologies, Inc. concerning your scooter.

Model _____ Serial Number _____

Date of Purchase _____ Body Color _____

Provider Name _____

Company _____

Address _____

City _____ State _____ Zip _____

Please remember to fill in and return your warranty registration card in the back of this manual.

IV. YOUR GOLDEN Buzzaround Luxury™

Specifications	GB119	GB149
Medicare Code	Not Coded	Not Coded
Weight Capacity	375 lbs.	375 lbs.
Drive Wheels	Rear	Rear
Maximum Speed	5 mph	5 mph
Operating Range ¹	18 miles	18 miles
Ground Clearance Under Deck	4 in.	4 in.
Turning Radius	40.5 in.	56 in.
Type Batteries	2-U1	2-U1
Colors	Red, Blue	Red, Blue
Freewheel Mode	Yes	Yes
Electronic Speed Control	Yes	Yes
Electro-Mechanical Brakes	Yes	Yes
Charger	Off-board	Off-board
Controller	90 Amp	90 Amp
Length	45 in.	44 in.
Width	22 in.	22 in.
Height (ground to top of seat back)	37-39 in.	37-39 in.
Height (ground to top of headrest)	43.5-48.5 in.	43.5-48.5 in.
Ground to Top of Seat	21.75-23.75 in.	21.75-23.75 in.
Ground to Top of Deck	6.75 in.	6.75 in.
Deck to Top of Seat	15.75-17.75 in.	15.75-17.75 in.
Number of Seat Height Positions	3	3
Size of Increments	1.0 in.	1.0 in.
Front Axle to Rear Axle	34.7 in.	33.8 in.
Weight of Unit: (assembled)	168 lbs.	174 lbs.
Front Half	50 lbs.	57 lbs.
Rear Half	36 lbs.	36 lbs.
Seat with Arms	32 lbs.	32 lbs.
Batteries	24 lbs./Pack x 2 = 48 lbs.	24 lbs./Pack x 2 = 48 lbs.
Tires:		
Standard Tires	Solid, low profile	Solid, low profile
Optional Tires	Not Available	Pneumatic ³
Front	9 in.	9 in.
Rear	9 in.	9 in.
Rear Anti-Tip Wheels	2 in.	2 in.

IV. YOUR GOLDEN Buzzaround Luxury™

Specifications	GB119	GB149
Standard Seat:	Highback	Highback
Seat Back Height (no headrest)	17 in.	17 in.
Seat Back Height (with headrest)	23-26 in.	23-26 in.
Width x Depth	17" x 16"	17" x 16"
Color	Black Vinyl	Black Vinyl
Optional Seat:	Stadium Style	Stadium Style
Seat Back Height (no headrest)	16 in.	16 in.
Seat Back Height (with headrest)	N/A	N/A
Width x Depth	20" x 18"	20" x 18"
Color	Black Vinyl	Black Vinyl
Warranty:		
Frame	Lifetime Warranty	Lifetime Warranty
Drive Train	2 Year Warranty	2 Year Warranty
Electronics ²	2 Year Warranty	2 Year Warranty
Batteries	1 Year Warranty	1 Year Warranty

Notes:

1. Battery range at 200 lbs., and will vary due to rider weight, drive surface, terrain, and battery type.
2. Electronics warranty excludes batteries.
3. Inflate tires to the maximum tire pressure recommended by the tire manufacturer.

Literature is current at the time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.

V. OPERATION



Figure 5. Buzzaround Luxury™ Delta Tiller Control Panel

CONTROL PANEL

Speed Adjustment Dial

- The speed adjustment dial allows you to preselect the maximum top speed of your Golden Buzzaround Luxury. See figure 5.
- Turn the dial counterclockwise to operate your scooter at the lowest speed level.
- Turn the dial clockwise to increase the operational speed of your scooter.
- We recommend that you select a speed setting at which you feel comfortable, safe, and in control of your Buzzaround Luxury.

Key Switch

This switch turns the power on and off to your Buzzaround Luxury. See figure 5.

- Insert the key into the key switch.
- Turn the key 90 degrees clockwise to turn on your scooter's power. The LEDs on the battery meter will light up indicating that power is on.
- Turn the key back to the vertical position to turn off your scooter's power.

CAUTION

Always make certain that the key is removed from the key switch before getting on or off your Golden Buzzaround Luxury or before lowering the tiller completely.

Horn Button

Press this button to sound the horn. See figure 5.

V. OPERATION

DC Battery Voltage Meter

This meter (see figure 5 on page 16) shows the amount of voltage in your Buzzaround Luxury's batteries.

- Green indicates a full charge on the batteries.
- Yellow indicates caution and approximately a one-half charge of the batteries.
- Red indicates that the batteries have less than one-half of a charge remaining.

See section X. "Care and Maintenance" for instruction concerning charging the batteries.

Delta Tiller Throttle Control Lever

The Delta Tiller enables you to use only one hand (either the right or left) to operate both the speed and the direction of the Buzzaround Luxury. This self-centering lever system controls your Buzzaround Luxury's speed (up to the maximum speed set by the speed adjustment dial) and its forward and reverse direction. See figure 5.

To drive forward use either of the following:

- Use your right hand fingers to pull back on the right side of the throttle control lever.
- Use your left thumb to push the left side of the throttle control lever.

CAUTION

Always bring your Buzzaround Luxury to a full stop before changing direction from forward to reverse, or from reverse to forward.

To drive in reverse use either of the following:

- Use your left hand fingers to pull back on the left side of the throttle control lever.
- Use your right thumb to push the right side of the throttle control lever.

When the throttle control lever is completely released, it automatically returns to the center "stop" position and engages your scooter's brakes, bringing you to a complete stop. You will hear a "click" when the park brake engages.

Control Options

Golden Technologies, Inc. has made provision for the individual needs and abilities of the operators of the Golden Buzzaround Luxury.

Your authorized Golden Technologies, Inc. provider can reverse the throttle control lever controls so that when you pull on the left side of the lever, your scooter will move in the forward direction, and when you pull on the right side of the lever your scooter will move in the reverse direction.

DELTA TILLER ADJUSTMENT

The tiller on your Buzzaround Luxury scooter employs an angle adjustment lever to provide you with multiple tiller angle settings. To set the tiller to your liking, follow these directions.

1. Squeeze the tiller angle adjustment lever. See figure 4 on page 14.
2. Position the tiller to a comfortable driving position.
3. Release the lever to lock the tiller.

V. OPERATION



Figure 6. 1 AMP USB Charging Port

DASH/LIGHTS
1 AMP FUSE
USB PORT
2 AMP FUSE
XLR PORT
10 AMP FUSE



Figure 7. Fuse Panel

USB Disclaimer: The supplied USB charging port is only intended to be used with appropriate smart devices such as phones or tablets. Please refer to device manufacturer's charging requirements before using. Golden Technologies excludes liability for any damages or data loss regarding the charging device on a computer or other electronic device connected to the USB port.



WARNING Fuses **MUST** be replaced with same type and rating.

SEATING



Figure 8. Armrest Width Adjustment



Figure 9. Armrest Angle Adjustment

Armrest Width Adjustment - See figure 8

To adjust the armrest width:

1. Loosen the adjustment knobs at the rear of the seat.
2. Slide the armrests in or out to the desired width.
3. Tighten the adjustment knobs.

Flip-up Armrests

Pull up on the end of either armrest to flip it up for easy transfer on and off your Golden Buzzaround Luxury™. See #1 on figure 9.

Armrest Angle Adjustment

Loosen the locknut by turning it counter-clockwise.

To increase the armrest angle, turn the armrest angle screw counterclockwise.

To decrease the armrest angle, turn the armrest angle screw clockwise.

See #2 on figure 9.



WARNING Be sure adjustments knobs are **securely tightened against the inserted armrest tube. Do not overextend the armrest.**

V. OPERATION



Figure 10: Seat Rotation Adjustment

Seat Rotation Lever

The seat can be rotated through 360° and locked at any 45° position.

1. Pull up on the seat rotation lever and rotate the seat to the desired position. See figure 10.
2. Release the handle to lock the seat at any 45° position.



Figure 11



Figure 12

Seat Height Adjustment

To adjust the seat height:

1. Fold the seat back down, grasp the seat on both sides and lift the seat up and off of the seat pedestal. See Figure 11. (*Note: Seat arms can be removed to reduce the seat weight*).
2. Remove both battery packs by lifting them up, one at a time. See figure 12.
3. Remove the nut and bolt that holds the seat post in the seat pedestal.
4. Position the seat post at the desired height while aligning the holes in the seat post and the seat pedestal.
5. Reinstall the nut and bolt.
6. Reinstall both battery packs.
7. Place the seat back onto the seat post and push down to lock it into place. Move the seat from side to side to make sure it is locked into place.

V. OPERATION



WARNING

Do not sit in or attempt to move your Buzzaround Luxury unless the frame and rear carriage are securely engaged. Accidental separation of the frame from the rear carriage may occur and cause injury or damage.



WARNING

Never sit on your Buzzaround Luxury in freewheel mode.



WARNING

Never disengage the park brake while you are operating your Buzzaround Luxury.



WARNING

Always make certain that your scooter park brake is engaged before inserting the key into the key switch and turning it to the on position.

FREEWHEEL MODE

Your Buzzaround Luxury™ is equipped with a freewheel lever that can set your scooter in or out of freewheel mode.



Figure 13: Freewheel Lever (Drive Position)

To disengage the park brake and put your scooter in freewheel mode:

Push forward on the freewheel lever towards the letter N (Neutral Position). See figure 13.

To engage the brake and take your scooter out of freewheel mode:

Pull back on the freewheel lever towards the letter D (Drive Position). See figure 13.

NOTICE

Pushing the scooter too quickly in freewheel will cause the motor to act as a generator causing the scooter to become difficult to push. Freewheel slowly.

VI. RIDING YOUR Buzzaround Luxury™

MOUNTING

CAUTION

For your first driving session, make certain that your Buzzaround Luxury starts out on a level surface and that you will continue driving on a level surface.

BEFORE GETTING ON YOUR BUZZAROUND LUXURY™

1. Check to be certain that the key has been removed from the key switch. See section V. "Operation" in this manual. This will eliminate the possibility of accidentally activating the paddle controls and causing injury to yourself and/or others.
2. Check to be certain that your Buzzaround Luxury freewheel lever is in the engaged position. See section V. "Operation" in this manual.
3. Flip up the armrest(s).



WARNING

When getting on or off your Buzzaround Luxury, keep your weight toward the middle of the deck. Putting most or all of your weight on the edge of the deck may cause an unstable condition.

GETTING ON YOUR BUZZAROUND LUXURY™

1. Position the seat for safe and easy mounting. Ensure that the seat is facing forward. See Section V. "Operation."
2. Return tiller to full upright position if necessary.
3. Carefully place one foot on the approximate center of the deck and seat yourself comfortably and securely on the seat.
4. Flip down the armrest(s).
5. Place the key into the key switch.
6. Rotate the key clockwise to the on position. See figures 14 and 15.



Figure 14. Key Switch (Off)



Figure 15. Key Switch (On)

VI. RIDING YOUR Buzzaround Luxury™

BASIC DRIVING

1. Make certain that you are seated safely and properly on your Golden Buzzaround Luxury.
2. Turn the speed control dial fully counterclockwise to its slowest setting.
3. Insert the key into the key switch—if you have not already done so.
4. Turn the key clockwise to the “On” position. See figure 15 on page 22.
5. Place your hands on the handgrips.
 - If you wish to drive forward, pull back on the right side of the throttle control lever (or push the left side of the throttle control lever forward.)
 - If you wish to drive in reverse, pull back the left side of the throttle control lever (or push the right side of the throttle control lever forward.)
6. Pull on the throttle control lever to gently accelerate your Buzzaround Luxury.
7. Release the throttle control lever to allow your Buzzaround Luxury to come gently to a full stop.
8. Practice these two basic functions until you feel that you have control of your scooter.

Steering

Steering your Golden Buzzaround Luxury is easy and logical.

1. With both hands on the handgrips of the tiller (see figure 4 on page 14), turn the tiller to the right to travel to the right.
2. With both hands on the handgrips of the tiller (see figure 4 on page 14), turn the tiller to the left to travel to the left.
3. Make certain to maintain sufficient clearance when turning your scooter so that the rear wheels will clear any obstacle.



WARNING

Turning your scooter too sharply at too high a speed may create a situation where one of the rear wheels will leave the ground. This may increase the possibility of tipping the scooter. Avoid this danger at all times by decelerating and steering a wide arc around corners and obstacles.

VI. RIDING YOUR Buzzaround Luxury™

Steering in a Tight Spot

If you must steer in a tight spot, such as entering a doorway or when turning around:

1. Bring your Buzzaround Luxury to a full stop.
2. Set speed at the lowest setting.
3. Turn the tiller to the direction in which you wish to drive.

For more information on dealing with obstacles and tight spots, see “Control Through Tight Spots” on the next page.

Steering in Reverse

Backing up your Buzzaround Luxury requires attention to what you are doing.

Use extreme caution when operating your scooter in reverse.

Plan your route to avoid getting into a difficult situation.

Make sure your proposed route is clear of people, pets, and obstacles.

1. Use your right fingers to push the throttle control lever or use your left fingers to pull back on the throttle control lever. See figure 5 on page 17.
2. Turn the tiller to the left to drive in reverse to the left.
3. Turn the tiller to the right to drive in reverse to the right.

NOTICE

Your Buzzaround Luxury’s speed in reverse is fifty percent of the speed set at the speed control dial.



WARNING

Prevent injury! Do not drive down an incline in reverse.

VI. RIDING YOUR Buzzaround Luxury™

CONTROL THROUGH TIGHT SPOTS

As you use your Golden Buzzaround Luxury to greatly increase your mobility, you will undoubtedly encounter some obstacles that will require practice to negotiate smoothly and safely. Below are some common obstacles that you may meet during the daily use of your scooter. Listed with those obstacles are some driving tips that should help you conquer those obstacles. Learn and follow those tips, and with surprising ease you will soon be in control of your Buzzaround Luxury as you maneuver it through doors, up and down ramps, up and over curbs, through grass and gravel, and up and down inclines.

Ramps

When proceeding up any ramp, curb, or incline:

- Lean forward in your seat to move your center of gravity forward for maximum stability and safety. If the ramp has a switchback, good cornering ability is required.
- Maneuver your Buzzaround Luxury so that the front wheels take wide swings around the corners of the ramp.
- Doing this will allow your scooter's back wheels to follow a wide arc around the corner, staying clear of obstacles.

If you must stop your Buzzaround Luxury while driving up a ramp:

- Starting up again simply requires that you apply gentle and steady forward power pressure to the throttle control lever.
- Accelerate gently after stopping on any incline.

Driving down a ramp:

- Keep your Buzzaround Luxury's speed control dial set fully counterclockwise at the slowest speed setting.
- If you must come to a stop, release the throttle control lever slowly and smoothly.



WARNING

Prevent injury! Do not drive down an incline in reverse.

Curbs

- See specifications for maximum curb height.
- Always use caution when negotiating any curb.
- Go up or down a curb head on at a direct 90° angle.
- Approach and negotiate the curb so that both back wheels of your scooter go over the curb at the same time.
- Never negotiate inclines or curbs by traversing them. Doing so may cause the scooter to tip over.
- Go down a curb slowly to avoid a jarring bump. Use as little power as possible.

Grass and Gravel

Your Golden Buzzaround Luxury performs admirably on grass, gravel, and hills, but you must follow the operational parameters presented in this manual. Refer to the "Safety" section. If you are unsure about any situation, avoid it. Common sense is your best protection.

- Feel free to use your Golden Buzzaround Luxury on lawns or in park areas.
- Avoid long or high grass, which may wrap around your scooter's axles.
- Avoid loose gravel.

VII. DISASSEMBLY

DISASSEMBLY

Your Buzzaround Luxury scooter is designed to be disassembled easily in order to be transported. When disassembling your scooter, please be sure to follow these simple instructions.

Remove the Seat

1. Place your scooter on a level surface.



WARNING

Remove the key from the key switch before you begin to disassemble your Buzzaround Luxury.

2. (Optional) Loosen the armrest adjustment knobs and remove the arms from the seat. This will reduce the weight of the seat making it easier to lift.



WARNING

Never pull on the arms or use them as a handle to move the scooter as they could inadvertently pull out of the receiver bar causing injury.

3. Grip the seat on opposite sides and, with a firm grip; pull the seat straight up. See figure 16.



Figure 16. Removing the Seat



Figure 17. Removing the Battery Pack

Remove the Battery Packs

4. Lift both battery packs off the scooter one at a time. See figure 17.

Remove the Front Basket

5. Remove the front basket by pulling up.

VII. DISASSEMBLY

Fold Down and Lock the Tiller

6. Push in and turn the tiller lock clockwise 90 degrees. See #1 on figure 18. This will lock the front wheel to keep it from turning side to side to help control while carrying.

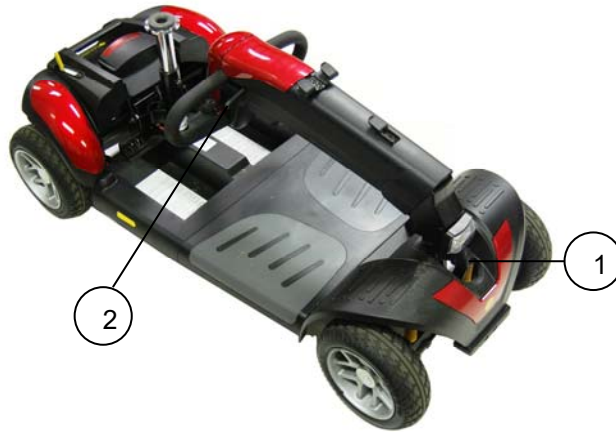


Figure 18. Lowered Tiller

7. Squeeze the tiller angle adjustment lever and fold the tiller down flat. See #2 on Figure 18. Then release the lever to secure the tiller in its new position.

Separate the Drivetrain from the Frame



WARNING! Be careful not to pinch your fingers between the frames. This area is a pinch-point and requires your full attention.

8. Pull up on the drivetrain release lever with one hand to unlock the front and rear frame sections. See Figure 19.

9. Pull the front frame up and off of the drivetrain. See figure 19.



Figure 19. Separating the Drivetrain and Frame

VIII. ASSEMBLY



WARNING! Be aware of any possible pinch-points when connecting the rear frame to the front frame. Please note warning labels at pinch-points during assembly.



WARNING! Pinching and crushing hazard! Increased hazards due to pinching or crushing result due to the component weight (such as batteries) during preparation for transport and maintenance work. Always carry out any work to be done with great care. Always try to get help from a second person, especially when stowing parts for transport.



WARNING! Pinch Point: A pinch point is an area where you **MUST** consider your fingers to be sure they **DO NOT** become pinched between two metal or hard plastic areas. The greatest chance of pinch points becomes possible during the assembly of the Buzzaround Luxury scooter.



WARNING

INJURY HAZARD DUE TO IMPROPER ASSEMBLY!

Ensure that all components in the scooter have been correctly assembled. During assembly, please check that all connecting/locking devices are holding in place.



Figure 20. Buzzaround Luxury™ Components (Model GB149)

VIII. ASSEMBLY

Reinstall the Drivetrain

1. Line up the frame and the drivetrain. See figure 21.



WARNING! WATCH FINGER PLACEMENT! BE CAREFUL NOT TO PINCH YOUR FINGERS!

2. Lower the frame onto the drivetrain until the drivetrain release lever engages. You should hear an audible click. See figure 22.

Install the Battery Packs

3. Gently lower the battery packs one at a time onto the frame. See figure 23.

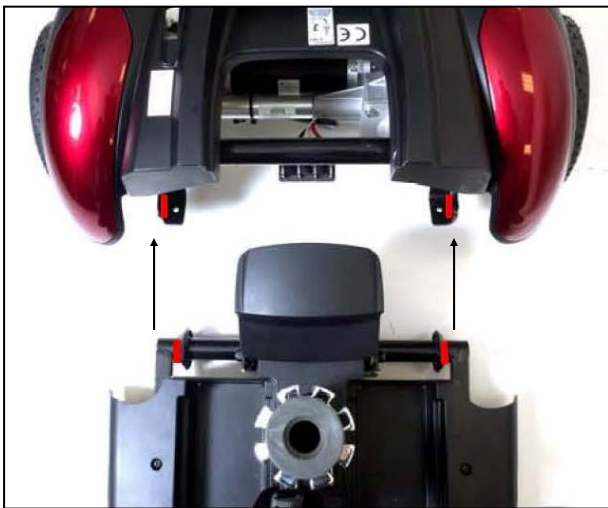


Figure 21. Drivetrain Alignment



Figure 22. Connecting Frame and Drivetrain



Figure 23. Battery Pack Installation

VIII. ASSEMBLY

Raise the Tiller

4. Squeeze the tiller angle adjustment lever and raise the tiller. See figure 24.
5. Release the adjustment lever to lock the tiller in the desired position.



Figure 24. Tiller Raised

Install the Front Basket

6. Place the basket onto the tiller bracket and push down to lock it in place. See figure 24.

Unlock the Tiller

7. Turn the tiller lock knob counterclockwise 90 degrees to the unlocked position. See figure 24.



WARNING IF YOU DON'T UNLOCK YOUR TILLER LOCK KNOB, YOUR TILLER WILL NOT TURN! THIS COULD CAUSE INJURY.

VIII. ASSEMBLY

Install the Seat

8. Hold the seat on an angle so you can see the seat swivel post underneath.
9. Place the seat swivel post into the seat post as shown in figure 25.
10. Push down on the seat to lock the seat in place. Move the seat from side to side to make sure it is locked in place.



Figure 25. Installing the Seat



WARNING! Pinch Point! Keep hands and clothing clear of the seat swivel post and seat post.



WARNING **Accident hazard due to non-engaged seat! Ensure that the seat is properly engaged by turning the seat slightly left and right.**

Accident hazard due to incorrect seat position! Always turn the seat to face forwards and engage it before driving. If the seat is in any other position, the possibility of tipping is increased. Before turning the seat, always ensure that the scooter is on an even and solid surface.

IX. TRANSPORTING YOUR Buzzaround Luxury™



Figure 26. Buzzaround Luxury™ (Model GB119)

NOTICE

All Golden Buzzaround Luxury scooters can be equipped with docking devices for loading onto a vehicle by means of a mechanical lift or hoist. Contact your Golden Technologies, Inc. provider for more information concerning docking devices and scooter lift devices.

To Transport Your Buzzaround Luxury:

- Disassemble or fold down the seat and tiller (see section VII. “Disassembly”) as close as possible to the loading area of the vehicle to be used for transporting your scooter.
- The degree to which you disassemble your Buzzaround Luxury or fold the seat and tiller depends on the size and shape of the cargo area of the transporting vehicle.
- Moving blankets or other forms of padding may be used to protect your Buzzaround Luxury during transportation.

CAUTION

- **Do not** lift your Buzzaround Luxury by its plastic body parts or dash. Breakage of these parts is not covered by the warranty.
- **Do not** use the tires or wheels to lift the rear carriage. The carriage may rotate and cause injury or damage.



WARNING

If your scooter and its components are not properly and securely stowed during transportation, the scooter or the components may move or become airborne and cause injury or damage.



WARNING

Never sit on your Buzzaround Luxury when it is being transported. Secure with an approved tie-down system.



WARNING

If transporting your scooter on an outside lift, remove seat or make sure the seat is fasten down securely with an approved tie-down system supplied by your lift manufacturer while transporting your scooter. Failure to do so could result in the seat coming loose and becoming airborne. This situation could cause problems for the traffic behind you and possibly result in an injury.

X. CARE AND MAINTENANCE

ROUTINE MAINTENANCE

Your scooter requires a minimal amount of care and maintenance. If you do not feel confident in your ability to perform the maintenance listed below, you can schedule an inspection and maintenance with your provider.

- ✓ Make sure your batteries are fully charged daily. See “Battery Charging”.
- ✓ Regularly inspect your scooter’s tires for signs of wear.
- ✓ Regularly inspect all exposed harnessing and connections for wear, damage, and corrosion. Have your provider repair or replace any damaged harness or connector.
- ✓ Periodically check all fasteners for tightness, wear, and corrosion. If any of these are present, turn off the scooter and contact your provider.



WARNING

Be careful of pinch points and sharp metal parts (wear gloves).

- ✓ Keep all the electrical components (charger, console panel, and electronics) free from moisture. If any of these items do become exposed to moisture, let them dry thoroughly before attempting to operate your scooter again.
- ✓ All bearings are permanently lubricated and sealed. No additional lubrication is required.

CLEANING YOUR Buzzaround Luxury™

Tires

- Clean the tires with ordinary kitchen-type cleaners and a damp cloth.
- Do not use solvents on the tires. Solvents can cause the tire material to break down or become too soft.
- **Optional Pneumatic tires only!** Check tire pressure every 3 weeks. For the proper tire pressure and inflation procedure, see page 36.

Body

The body (rear fender shrouds (2), nose piece (1), and the console shroud) of your scooter is formed ABS plastic and is painted for style and durability. Each Buzzaround Luxury comes standard with two sets of quick snap color shrouds in Red and Blue. Colors can be changed by gently snapping off the pieces.

NOTICE

Do not use excess water or pressure when cleaning your scooter.

- **Never** take your scooter through a car wash.
- **Never** use a high pressure hose to wash your scooter.
- Use cool water mixed with a mild soap to remove dirt and oils.
- Dry with a clean soft cloth.
- Buff by hand with a soft cloth.
- **Do not** use a sharp object when removing shroud pieces.
- Clean the body of your scooter with a damp cloth.

Vinyl Seat

Clean with a mild soap or mild detergent and a damp cloth. A vinyl cleaner may also be used.

X. CARE AND MAINTENANCE

IMPORTANT! New batteries **MUST** be fully charged prior to the initial use of your scooter. **Charge new batteries for 12 hours even if the battery gauge already shows a full charge. This is to condition your new batteries for maximum performance.**



WARNING! Corrosive chemicals contained inside the battery.
Explosive battery conditions exist.

WARNING! Charger surface can become HOT.
Electrical hazard!

BATTERIES AND CHARGING

Battery maintenance is the most important part of maintaining your scooter. Keeping the batteries fully charged helps to extend battery life. Use the following guidelines to help keep your batteries in optimum condition.

- For daily use, keep batteries fully charged. We recommend that you plug in the off-board charger after each use and charge 6-8 hours.

“For best results please charge your batteries overnight after every use.”

- If you are not going to use the scooter for more than a week, **fully charge** the batteries and then remove the battery packs from the scooter. See figure 17 on page 26.

Charging Guidelines Checklist to Maximize Battery Life

- ✓ Use only the automatic off-board charger supplied for all routine charging.
- ✓ Never use an automotive or wet type charger.
- ✓ Avoid deep discharges and never drain the batteries completely.
- ✓ Do not leave batteries in a low state of charge for extended periods.

Charge a discharged battery as soon as possible.

- ✓ Fully recharge batteries regularly.
- ✓ Always store batteries fully charged.
- ✓ Check stored batteries once a month and recharge as necessary.

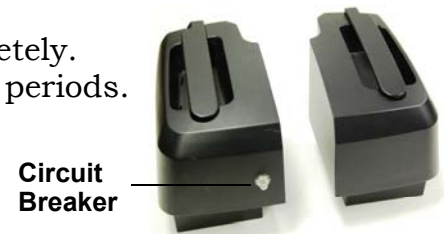


Figure 27. Battery Pack

To Charge the Batteries:

The charger for all Buzzaround Luxury models is an off-board charger. To charge the batteries, plug the charger into the charger port at the center of the tiller. See figure 28. When the batteries are fully charged, disconnect the charger power cord from the wall outlet and from the charging socket.

Note: Your scooter will not move while the charger remains connected.

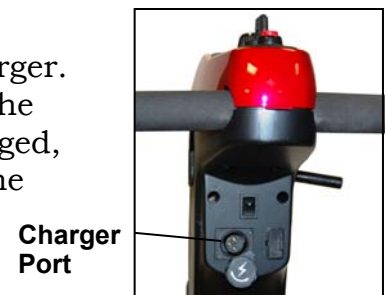


Figure 28. Tiller

Battery Maintenance

MANDATORY! Use only AGM or GEL-cell SLA deep-cycle batteries for operating the Golden Buzzaround Luxury scooters.

- These batteries are maintenance free.
- There is no danger of spillage or leakage, so these batteries are safely transportable on aircraft, buses, trains, etc.
- By following the procedures set out in this manual, you can expect extended life from your Golden Buzzaround Luxury's batteries.

X. CARE AND MAINTENANCE

IF YOUR Buzzaround Luxury™ DOES NOT OPERATE

- Make certain the freewheel lever is set to the engaged position. See page 21.
- Check the main circuit breaker. If necessary, reset the circuit breaker. See figure 27 on page 34.
- Make certain that the speed adjustment dial is at the desired setting.
- Make certain that the key switch is set to the “On” position.

If none of the above procedures solve the problem, contact your authorized Golden Technologies, Inc. provider.

MAIN CIRCUIT BREAKER

The main circuit breaker reset button is located on the battery pack. See figure 27 on page 34.



WARNING! Do not attempt electrical repairs.

Consult Golden Technologies, Inc. or your authorized Golden Technologies, Inc. provider. If, for no apparent reason, your Buzzaround Luxury stops operating; the main circuit breaker may have tripped.

Possible causes for the main circuit breaker to trip:

- Driving up a long, steep hill.
- Driving over a curb.
- Run-down batteries.
- Exceeding the maximum weight capacity of the scooter.

As you operate your Buzzaround Luxury, battery voltages go down and battery current must rise to satisfy the demands of the motor and of other electrical devices operating on your scooter. This can cause a heavy current draw that will trip the main circuit breaker.

Remedy:

- Recharge your Buzzaround Luxury's batteries. See “Charger and Batteries” in this section of this manual.
- If the problem continues, have your authorized Golden Technologies, Inc. provider load test the batteries.
- If the batteries are good, the charger may be the problem. Consult your authorized Golden Technologies, Inc. provider.

Resetting the Main Circuit Breaker

1. If the main circuit breaker trips as a result of run-down batteries or because of a temporary overload, reset the circuit breaker.
2. Wait ten minutes or so for the motor control board to return to the normal operating temperature range.
3. Make certain that the key switch is set to the “Off” position.
4. Press in on the main circuit breaker reset button. See figure 27 on page 34.

NOTICE

If the main circuit breaker continues to trip, there is probably an underlying electrical fault that needs attention.

- Contact your authorized Golden Technologies, Inc. provider.
- **Do not** keep resetting the main circuit breaker without correcting the underlying electrical problem.

X. CARE AND MAINTENANCE

PROPER TIRE INFLATION - OPTIONAL PNEUMATIC TIRES ONLY!



WARNING Golden Technologies, Inc. does not recommend that you disassemble the tire from the split rim wheels. Consult your authorized Golden provider for assistance with this matter. If you must disassemble a pneumatic tire from the split rim, the tire must be completely deflated. If the tire is not completely deflated, serious injury may occur.

NOTICE

Under-inflated tires will cause significant loss of operational range and will also cause the motor control board to overheat.

To inflate a tire:

1. Remove the valve stem cap by turning it counter-clockwise.
2. Inflate the tire to the maximum tire pressure recommended by the tire manufacturer.
3. Replace and turn the valve stem cap clockwise to secure it.
4. Repeat steps 1-3 for each tire.



WARNING Do not exceed the maximum tire pressure recommended by the tire manufacturer .

DIAGNOSTICS

The diagnostics feature of your Golden Buzzaround Luxury microprocessor based motor control board functions as follows:

- Any fault condition with the controller or with an associated system will cause the battery meter to flash and the horn to beep.
- The flashing of the meter and beeping of the horn occurs as one flash/beep or a series of flashes/beeps separated by a two-second pause.
- The number of flashes/beeps in each series is referred to as the “Flash/Beep Code.”
- The flash/beep code indicates the nature of the condition or fault.

Flash/Beep Codes

1 Flash/Beep

- Indicates that the battery voltage (with the throttle control lever in the neutral position) is below the required operating voltage of the scooter.
- The batteries need to be charged.
- Your Buzzaround Luxury will continue to operate at reduced speeds.
- As the voltage drops you will experience a power loss.

Remedy: Charge your Buzzaround Luxury’s batteries. See “Batteries and Charging” in this section.

2 Flashes/Beeps

- Indicates that the battery voltage is too low for the scooter to operate.
- At this voltage your scooter will cease operation.
- If you have charged the batteries and the condition continues, one or both of the scooter’s batteries may be at fault.
- The continuance of this condition after you have charged your Buzzaround Luxury’s batteries may also indicate a problem with the battery charger.

X. CARE AND MAINTENANCE

3 Flashes/Beeps

Consult your authorized Golden Technologies, Inc. provider.

4 Flashes/Beeps

- Indicates that the motor control board is overloaded and overheated.
- An overload can occur if you have been driving your Buzzaround Luxury for an extended period of time up an incline above the recommended 6 degrees.
- An overload can occur if your Buzzaround Luxury is carrying a payload that is higher than the recommended weight capacity of 375 pounds.
- Your scooter will not operate until the motor control board cools back to its operational temperature range.
- ***Optional Pneumatic Tires only!*** Check to make sure each tire is properly inflated. See "Proper Tire Inflation" on page 36.

Remedy:

- Turn off your Buzzaround Luxury's key switch and allow the scooter to sit for 10 minutes.
- Consult your authorized Golden Technologies, Inc. provider.

5 Flashes/Beeps

- Indicates a brake problem.
- The freewheel lever may be in the disengaged position. See page 21.
- The brake or the brake wiring may be damaged.
- For safety reasons, your scooter was designed to cease operation until the problem is corrected.

Remedy:

- Make certain that the freewheel lever is in the engaged position. See page 21.
- Turn the key switch to the "Off" position to stop the beep code.
- Set the freewheel lever to the "engaged" position.
- Turn the key back to the "On" position.
- If the above remedies do not solve the problem, contact your authorized Golden Technologies, Inc. provider.

6 Flashes/Beeps

- Indicates that the throttle control lever was not in the neutral position when the key switch was turned to the "On" position.
- May indicate that the throttle control lever is out of adjustment.

Remedy:

- Turn the key switch to the "Off" position.
- Make certain that the throttle control lever is in the neutral position.
- Turn the key switch to the "On" position.
- If the above procedure does not solve the problem, contact your authorized Golden Technologies, Inc. provider.

X. CARE AND MAINTENANCE

7 Flashes/Beeps

- Indicates a problem with the throttle control lever.
- May also indicate a problem with the potentiometer (speed control.)
- Your scooter will not operate until the problem is resolved.

Remedy:

- Make certain that all electrical connections are firmly and correctly joined.
- If the connections are all firmly joined and the problem is not solved, contact your authorized Golden Technologies, Inc. provider.

8 Flashes/Beeps

- Indicates a motor voltage problem.
- Your scooter will not operate until the problem is corrected.

Remedy: Contact your authorized Golden Technologies, Inc. provider.

9 Flashes/Beeps

- Indicates other internal errors or faults
- Your scooter will not operate until the problem is corrected


Remedy: Contact your authorized Golden Technologies, Inc. provider.

NOTICE


FOR YOUR CONVENIENCE, THESE DIAGNOSTIC FLASH/BEEP CODES ALONG WITH THE SERIAL NUMBER OF YOUR SCOOTER CAN BE FOUND UNDER THE BATTERY PACKS. THIS WILL ALLOW YOU TO ACCESS THIS IMPORTANT SERVICE INFORMATION AT ANY TIME. IF POSSIBLE, PLEASE REFER TO THIS INFORMATION BEFORE CONTACTING YOUR AUTHORIZED GOLDEN TECHNOLOGIES, INC. PROVIDER.

Non-Flash/Beep Code Faults

If the following faults occur, please check the fuses in the fuse panel located at the center of the tiller. Refer to page 19 for fuse replacement and warnings.

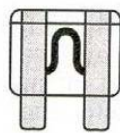

Good Fuse

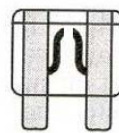

Bad Fuse

Typical Blade Style Fuse shown

Faults

- The dash does not illuminate when the key is turned to the “on” position. The scooter does not operate and the headlight/taillight do not illuminate.

Remedy: Check the top fuse (1 AMP) in the fuse panel. Replace the fuse with the same type and rating if the fuse is bad.

- The USB charging port will not charge your smart devices such as a phone or tablet.

Remedy: Check the middle fuse (2 AMP) in the fuse panel. Replace the fuse with the same type and rating if the fuse is bad.

- The XLR charging port is not functioning, causing the batteries not to charge when the off-board battery charger is connected.

Remedy: Check the bottom fuse (10 AMP) in the fuse panel. Replace the fuse with the same type and rating if the fuse is bad.

XI. WARRANTY



Consumer Limited Warranty *Scooters*

Golden Technologies has provided information throughout the owner's manual in regards to proper usage and care of your scooter. Failure to adhere to this information will void the warranty of the scooter, it's parts and the batteries supplied by Golden Technologies.

The following warranty is granted only to the INITIAL Consumer who has purchased our product and commences on the Date of Purchase by the Consumer from an Authorized Golden Technologies' Provider. The Consumer must have a warranty registration card filed with Golden Technologies to receive service.

The warranty covers the following scooter models unless otherwise noted:

- | | | |
|---|--|--|
| ◆ BUZZaround LT:
GB107 | ◆ Buzzaround Extreme:
GB118, GB148 | ◆ All Companion Models:
GC240, GC340, GC440 |
| ◆ BUZZaround XL Series:
GB117D, GB117S, GB117H, GB117Z
GB147D, GB147S, GB147H, GB147Z | ◆ BUZZaround Luxury:
GB119, GB149 | ◆ Avenger: GA541 |
| | ◆ LiteRider Series:
GL109, GL111, GL141 | ◆ Patriot: GR575 |

Lifetime Limited Warranty:

For the life of the scooter, from the date of purchase from an authorized Golden Technologies' dealer, in the event of defective materials or workmanship, Golden Technologies will repair or replace at our option any of the following structural frame components found to be defective by an authorized Golden Technologies representative:

- | | | |
|---------------|----------------|-------------|
| ◆ Steel Frame | ◆ Tiller Frame | ◆ Seat Post |
| ◆ Frame Welds | ◆ Front Fork | |

****The GR575 is warranted for FIVE YEARS***

Two-Year Warranty:

For the period of two years from the date of purchase, in the event of defective materials or workmanship, Golden Technologies will repair or replace at our option with any new or reconditioned components any of the following drive train components found to be defective by an authorized Golden Technologies representative:

- | | |
|-------------|---|
| ◆ Transaxle | ◆ Motor/Brake Assembly (electrical function only) |
|-------------|---|

Note: An increase in operating noise of the transaxle DOES NOT constitute a defect. With normal wear and tear, operating noise is expected to increase.

For the period of two years from the date of purchase, in the event of defective materials or workmanship, Golden Technologies will repair or replace at our option any of the following electronic assemblies found to be defective by an authorized Golden Technologies representative:

- | | | |
|-------------------------|--------------------------|--------------------|
| ◆ Electronic Controller | ◆ Dash Assembly | ◆ Wire Harness(es) |
| ◆ Battery Charger | ◆ Potentiometer Assembly | |

Note: While charging batteries, it is normal for the battery charger to heat up. Heat coming from the battery charger does not constitute a defect.

One-Year Limited Warranty

For the period of one year from the date of purchase, in the event of defective materials or workmanship, Golden Technologies will replace the **batteries supplied by Golden Technologies** if found to be defective by an authorized Golden Technologies representative.

XI. WARRANTY



Consumer Limited Warranty Scooters (continued)

One-Year Limited Warranty:

For the period of one year from the date of purchase, in the event of defective materials or workmanship, Golden Technologies will repair or replace at our option any of the following parts found to be defective by an authorized Golden Technologies representative:

- | | | |
|-------------------------|---------------------|----------------------|
| ◆ Bearings | ◆ Rubber Components | ◆ Plastic Components |
| ◆ Bushings | (excludes tires) | (excludes body) |
| ◆ Seat Swivel Mechanism | ◆ Armrests | |

Items not listed are covered at the discretion of Golden Technologies. Any accessories, standard or optional, supplied by Golden Technologies, are covered for a period of one year from the date of purchase with the particular unit.

Warranty Exclusions:

Golden Technologies does not provide warranty on any of the following items which may require replacement due to the normal wear and tear of day to day usage:

- | | | |
|-----------------------|-----------------|--------------|
| ◆ Tires and Tubes | ◆ Motor Brushes | ◆ Brake Pads |
| ◆ ABS Plastic Shrouds | ◆ Bulbs/Fuses | ◆ Upholstery |
| ◆ Armrest Pads | | |

This warranty also excludes the following:

- ◆ Loss or theft of components (ie: keys, baskets, etc.)
- ◆ Damage caused by:
 - Battery fluid spillage or leakage
 - Abuse, misuse, accident or negligence.
 - Improper operation, maintenance or storage.
 - Commercial use or use other than normal.
 - Repairs and/or modifications made to any part without the specific consent of Golden Technologies.
 - Exceeding the specified weight capacity of the unit.
 - Accessories other than those supplied or approved by Golden Technologies.
 - Failure to adhere to the product instructions
 - Acts of Nature, such as lightning strikes, earthquakes, etc.....
 - Circumstances beyond the control of Golden Technologies.

◆ ANY PARTS ALTERED OR REPAIRED BY UNAUTHORIZED PERSON.

Warranty Service:

Warranty Service must be performed by an authorized Golden Technologies representative. Golden Technologies reserves the right to replace warranted part(s) with refurbished or new part(s) at our discretion. All labor charges, service calls, transportation costs or any other charge(s) associated with the installation of any warranted part(s) are the responsibility of the consumer. Warranty is for the replacement of the parts only and does not include freight for the replacement parts. Consumers are not to return any item(s) to Golden Technologies without prior written authorization. Any damages incurred while warranted part(s) are in transport are the sole responsibility of the consumer.

There is no other express warranty.

Implied warranties, including those of merchantability and fitness for a particular purpose are excluded. Liabilities for consequential damages are excluded. This warranty gives you specific rights and you may also have other rights which may vary from state to state.

GOLDEN PRODUCT REGISTRATION

Register Now To Lock In Your Warranty and Receive Important Product Updates

SAVE TIME - REGISTER ONLINE: WWW.GOLDENTECH.COM

SERIAL NUMBER

**Most Lift Chairs Have On Controller Wire or Back Rail of the Metal Frame
Most Scooters & Powerchairs On Seat Post, Under/On Battery, or Under Carpet**

[illegible]**PRODUCT MODEL NUMBER**

Most Scooters with G (GA, GB, GC, GL, GR, etc) Most Powerchairs with GP

[illegible]

PURCHASE DATE

FIRST NAME

[illegible]

LAST NAME

[illegible]

ADDRESS

[illegible]

ADDRESS LINE 2

[illegible]

CITY

[illegible]

STATE / PROVINCE ZIP / POSTAL CODE

[illegible]

E-MAIL ADDRESS

[illegible]**PROVIDER NAME**[illegible]**PROVIDER'S ADDRESS**[illegible]

ADDRESS LINE 2

[illegible]

CITY

[illegible]

STATE / PROVINCE

[illegible]

ZIP / POSTAL CODE

[illegible]

COUNTRY

[illegible]

Product comments:

Product comments:
Rate Provider Experience (Circle One - 1 Bad, 5 Great)
1 2 3 4 5

Fold Here First



Place
Stamp
Here

Golden Technologies
401 Bridge St.
Old Forge, PA 18518



Fold Here Second



GB119



GB149



401 Bridge Street
Old Forge, PA 18518
Tel: 800-624-6374
Fax: 800-628-5165
www.goldentech.com