

L U X U R Y A D J U S T A B L E B E D



Shown with Premium Mattress

Sit Up



Lounge



Relax



Lay Flat



Rest

- Increase breathing capacity
- Improve circulation
- Relieve back pressure
- Minimize pressure points



3101 The ultimate sleep solution.



Rise. Roll. Rest. Ride.



YOUR LIFE IN MOTION

Bed Base Specifications

Massage System:	No
Wood Skirting:	No
Headboard Brackets:	1 Set Standard
Footboard Brackets:	Option
Half Length Rails:	Option
Bed Risers 2"-6":	Option
Wireless Remote System:	N/A
Weight Capacity (twin):	400 lbs
Weight Capacity (full):	600 lbs
Weight Capacity (queen):	600 lbs

Mattress Specifications



Premium:	
Quilt Pattern:	12" on Center
Ultra Quilt (layer 1):	1½" Ultra Quilt
Foam Support (layer 2):	1½" Visco Foam (Memory)
Foam Support (layer 3):	¾" High Density Foam
Posture Center (layer 4):	Yes
Edge Guards:	Yes
Cover Pattern:	Ivory

Size Specifications

Extra Long Twin Size	38" x 80"
Extra Long Twin Coil Count	234
Extra Long Full Size	54" x 80"
Extra Long Full Coil Count	338
Queen Size	60" x 80"
Queen Coil Count	390
Dual King* Size	76" x 80"
Dual King* Coil Count	568

Warranty

Bed Base	15 years
Premium Mattress	10 years

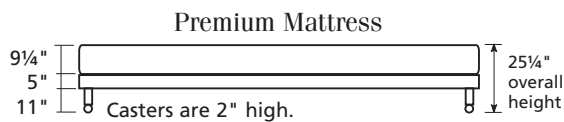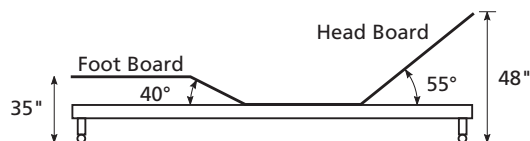

HAND CONTROL

Offering an exceptional value, the Golden Technologies standard model adjustable bed features independent head and foot operation with a state of the art hand control. Superior engineering and strong metal mechanisms offer smooth operation while a whisper quiet, low voltage motor powers the independent head and foot positioning.


- Low voltage electric motor system
- One master control wand for easy operation
- Available extra long twin, full, queen, and dual king
- Made in the U.S.A.

Solid Construction

Golden Technologies adjustable beds are manufactured at our factory in Old Forge, PA. Experienced technicians craft each adjustable bed by hand and our qualified quality control team inspects and approves each bed before shipping. Experience the difference with an adjustable bed from Golden Technologies.

Premium Mattress Construction/Selection

Our Premium, no-flip quilted mattress features 1½" Ultra Quilt topping with 1½" Visco foam, ¾" high density foam, posture center, edge guards, and a soft Ivory ticking.

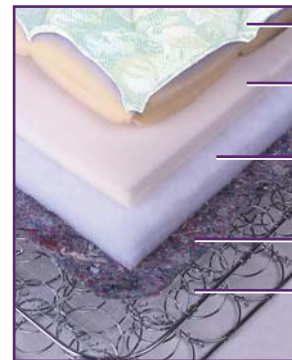

Shown with optional half rails.



Quality Engineering

Low voltage, quality motor system operates both the head and foot. Completely enclosed, single motor system design offering safe, durable operation. Dual motor massage system applied directly to the baseboard with a sound dampening application. Massage units are great for soothing the mind and body.

MADE IN USA

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.

* Dual King consists of 2 twin mattresses



Golden Technologies • 401 Bridge Street, Old Forge PA 18518 • Tel: (800) 624-6374 • Fax: (800) 628-5165 • www.goldentech.com