

M A X I C O M F O R T ® S E R I E S



PR761 MaxiComfort  
with Twilight Positioning  
shown in Shiitake

Introducing the  
Next Generation of  
Wellness Positioning

**NEW! Power  
Lumbar & Power  
Headrest  
STANDARD!**



Bourbon



Smoke



Brisa Coffee Bean\*



Brisa Shiitake\*



Brisa Iron\*

\* upcharge applies

**Pocket Coil Spring Seat Standard!**



*EZ<sup>™</sup>Sleeper*



**TWILIGHT**

**EZ Sleeper  
PR761-MLA**

The new EZ Sleeper with Twilight positioning technology features a modern, sleek design with a stylish track arm for ultimate comfort! Elevate and circulate in the PR761, the world's most advanced power positioning lift recliner. Golden's patented 5-motor system provides unique, rejuvenating positions along with power lumbar support and power headrest. Healthy benefits include being able to elevate the feet above the heart while experiencing a feeling of weightlessness. Twilight positions relieve pressure points and encourage proper posture by supporting the natural "S" shape of the spine.



Model Number	PR761-MLA
Size	Medium/Large
Weight Capacity	375 LBS
Positioning	Infinite
Twilight Technology	Yes
Chaise Seat	Yes
Pocketed Coil Seat	Yes
Back Type	Seam
Removable Back Option	No
Head Rest Cover	No
Arm Rest Covers	No
Overall Width	34"
Overall Height	46"
Width Between Arms	22"
Floor to Top of Seat	21"
Seat Depth	21"
Seat to Top of Back	30"
Seat to Top of Arm	8"
Distance Required from Wall Reclined	24"
Extension Reclined from Back of Seat to End of Foot Rest	43"
Life Time Warranty, Chair Frame, Lift Frame, Recline Mechanism	Yes
Warranty on Electrical parts	3 YR
Pro Rated Years 4 Thru 7	Yes
Warranty on Mechanical Labor	3 YR
Warranty on Electrical Labor	1 YR
Customizable *	No
Hand Control Routing	Inside Arm
Extra Pocket Option	Standard
Heat & Massage Option	No
Power Pillow Option	Standard
Power Adjustable Lumbar Option	Standard
AutoFlex Option	No
Foot Rest Extension Option	No
Straight Lift Option	No
USB Charging Port	Yes
Place of Origin	USA

\* See Customize Chart for Dimension limits  
 Sizes may vary depending on fabric, filling material, upholstery and carpet thickness. All measurements taken on concrete floor using levels and metal rulers.



### Smart Tek

Our luxurious lift and recline chairs are equipped with electronics to help assure a secure lift every time you need it.

### Battery Backup System

Batteries not included. Battery Backup System is designed to bring the chair from the reclined to the seated position in the event of a power outage.



Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.




Complies with  
 UL 73 and  
 CSA C22.2 No. 68



  
**Twilight**  
*ZZZ*  
**EZ Sleeper**  
**The World's Most Advanced Power Positioning Lift Recliner!**

**Our Patented Twilight Technology Now in a Unique 5-Motor System to Achieve Optimum Support**

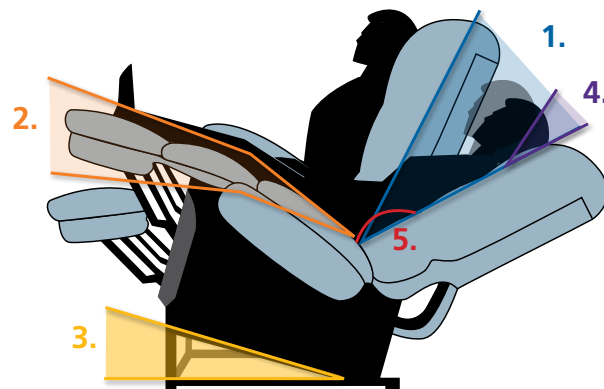

**Motor #1**  
 Adjusts back positioning

**Motor #2**  
 Controls the footrest movement and the patented movement of the rising seat box

**Motor #3**  
 Activates Twilight technology, a tilting movement that provides a weightless experience

**NEW! Motor #4**  
 Adjusts head support and placement.

**NEW! Motor #5**  
 Provides customizable lumbar support.

### HEALTHY BENEFITS:

- Elevates Feet Above Heart
- Supports Natural 'S' Shape of Spine
- Alleviates Pressure On Lungs
- Stylish New Track Arm Design
- Distributes Weight
- Available in Luxurious Brisa Fabrics With Breathable Technology
- Relieves Pressure Points
- Encourages Proper Posture