



GB148 Shown

Standard Colors:



Velvet Red
(Installed)



Coastal Blue
(Included in Box)

NEW Sporty Optional Colors:




Crisp White



Canary Yellow




Stone Grey

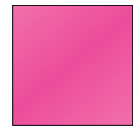

Orchid Pink



Army Green
(Matte)

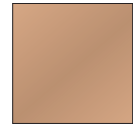

Desert Tan



Sporty New LED Running Light on the 4-W EX

Buzzaround 3-wheel GB118 4-wheel GB148



The Golden Buzzaround Extreme Scooter features the portability and easy disassembly of a lightweight scooter with the amazing benefits of a full-size scooter! Standard comfort features on the Buzzaround EX include a stylish quilted high back seat, an infinite adjustable tiller at your fingertips, and a USB charging port right under the dashboard. The 4-Wheel model features a sporty new LED V-shaped running light on the front for maximum safety and visibility. The Buzzaround EX also comes standard with 2-U1 batteries - industry leading and the largest battery in this category - for maximum operating range up to 18 miles. With a 350 lb. weight capacity, comfort spring suspension, and 40.5" turning radius, the Buzzaround EX offers a smooth ride. Go the extra mile on a BuzzAround Extreme Scooter!

Standard Features:

- 350 lb. weight capacity
- 4" ground clearance
- Full size & portable
- Larger 90 amp controller
- Larger 320 watt motor
- Larger 9" tires
- U1 battery for maximum range
- Ultrabright LED headlight & tail light
- Comfort suspension springs
- Full size, yet portable & easy to disassemble





Wireless disassembly in 4 easy steps!

Model:	GB118	GB148
Description	EX	EX
Class I Medical Device	Yes	Yes
Weight Capacity	350 lbs.	350 lbs.
Drive Wheels	Rear	Rear
Maximum Speed	5 mph	5 mph
Operating Range up to ¹	18 miles	18 miles
Under Deck Clearance	4"	4"
Turning Radius	40.5"	54"
Type Batteries	2-U1	2-U1
Colors	Red, Blue	Red, Blue
Freewheel Mode	Yes	Yes
Electronic Speed Control	Yes	Yes
Electro-Mechanical Brakes	Yes	Yes
Charger	Off-board	Off-board
Controller	90 Amp Dynamic R-Series	90 Amp Dynamic R-Series

Dimensions:		
Length	45"	44"
Width	22"	22"
Height (ground to top of back of seat)	35"-37"	35"-37"
Ground to Top of Seat	22.75"-24.75"	22.75"-24.75"
Ground to Top of Deck	6.75"	6.75"
Deck to Top of Seat	16"-18"	16"-18"
Number of Seat Height Positions	3	3
Size of Increments	1"	1"
Front Axle to Rear Axle	34.7"	33.8"

Weight of Unit: (assembled)		
Front Half	154 lbs.	161 lbs.
Rear Half	46 lbs.	53 lbs.
Seat w/Arms	36 lbs.	36 lbs.
Batteries	23 lbs.	23 lbs.
	24 lbs. per pack	24 lbs. per pack

Tires:		
Standard Tires	Solid, low profile	Solid, low profile
Optional Tires	N/A	N/A
Front	9"	9"
Rear	9"	9"
Rear Anti-Tip Wheels	2"	2"

Standard Seat:	Highback	Highback
Seat Back Height - no headrest	18"	18"
Width x Depth	18" x 16"	18" x 16"
Color	Black Vinyl	Black Vinyl
Optional Seat:	Stadium Style	Stadium Style
Back Height - no headrest	16"	16"
Width x Depth	20" x 18"	20" x 18"
Color	Black Vinyl	Black Vinyl

Warranty:		
Frame	Lifetime Warranty	Lifetime Warranty
Drive Train	2 Year Warranty	2 Year Warranty
Electronics ²	2 Year Warranty	2 Year Warranty
Batteries	1 Year Warranty	1 Year Warranty

¹ Battery range at 200 lbs. and will vary due to rider weight, drive surface, terrain, and battery type.
² Electronics warranty excludes batteries.



Step 1: Pop off the seat!



Step 2: Pull out the battery pack!



Step 3: Fold down tiller and lock in place!



Step 4: Pull forward on the drivetrain release lever to separate the front and rear frame sections!



Disassembles into 6 pieces!

4 wheel model shown



A battery charging port is conveniently located on the tiller under the control panel.



Phone: 800-265-9615
www.goldentechnologies.ca



www.facebook.com/goldentechcanada