



Standard Colors



Grey



Teal

**Remote
Control App
for Enhanced
Convenience**

**Folds Easily
Drives Remotely
Tight Turning
Radius**

Standard Features:

- Smart Control Mobile App
- Memory Foam Seat Cushion
- Padded Armrests
- Adjustable Recline Back
- Spacious Leg & Foot Area
- Length Adjustable Footrest
- Front Suspension
- Omni-Directional Wheels
- Airline Friendly, Lithium-Ion Battery
- Easy-to-Reach Storage Compartment
- System Monitoring Voice Alerts
- USB Charging Ports
- Front & Rear Lighting
- Water Resistant Backpack

**Exclusive
Omni-Directional
Wheels with
Self-Cleaning
V-Groove!**

**AllyPro
GP304**



Visit App Store
for Remote
Control



The Golden Ally Pro™ is next in line for our smart power wheelchair series. It features a modern design that applies robotics technology for easy folding and remote control. Upgrades from the Ally™ include 10" omni-directional front wheels allowing for more leg and foot room, the ability to climb larger obstacles, and rolling forward and sideways enabling the chair to turn in a tight radius. It offers great ride quality and capability, compliments of its large, pneumatic rear drive wheels and shock-absorbing front suspension. The two-position adjustable back angle provides long-lasting comfort. Arms are adjustable for additional seat width and feature padded armrests. Its comfortable seating and foot positions will also be a welcome change from the compromises of the typical power chair. The Ally Pro is perfect for small spaces indoors or open spaces outside. With a long-lasting lithium-ion battery, user-friendly interface, and Smart Control Mobile App, you can navigate any environment with ease and confidence. Move effortlessly with the Ally Pro, your perfect travel partner!




Model Number	GP304
Class I Medical Device	Yes
Weight Capacity	330 lbs.
Drive Wheels	Rear
Maximum Speed	5 mph
Operating Range ¹	13 miles
Ground Clearance	2.5"
Turning Radius	37"
Type Batteries	12 AH Lithium-Ion
Colors	Teal, Grey
Freewheel Mode	Yes
Electronic Speed Control	Yes
Electro-Mechanical Brakes	Yes
Charger	Off Board
Controller	40 Amp Remote
Dimensions:	
Length	39"
Width	25"
Width Between Arms	20" max.
Height (ground to top of back of seat)	37"
Ground to Top of Seat	20.5"
Weight of Unit: (assembled)	75 lbs.
Heaviest Piece	69.6 lbs.
Battery	5 lbs.
Tires:	
Drive Wheels	10" Pneumatic
Casters	10" Omni-Directional
Anti-Tip Wheels	2"
Standard Seat:	
Back Height - no headrest	19"
Width x Depth	17" x 16"
Color	Black Fabric
Footrest:	
Height adjustable	Yes
Warranty:	
Frame	5 Years
Drive Train	13 Months
Electronics ²	13 Months
Batteries	12 Months

¹ Battery range tested at 200 lbs. and will vary due to rider weight, drive surface, terrain and battery type.

² Electronics warranty excludes batteries.



Padded Armrests

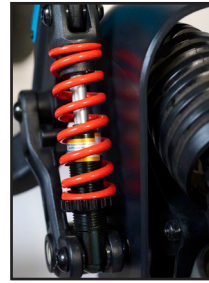

Front Suspension



Highly Visible Rear Lighting



Zippered Water Resistant Backpack



Omni-Directional Wheels with Self-Cleaning V-Groove



Pre-wired Controller for right or left side placement

Optional Accessories:



Extendable Headrest



Cup/Phone Holder



Heated Seat Control



Wireless Remote



Phone: 800-265-9615
www.goldentechnologies.ca



www.facebook.com/goldentechcanada